

head

head



head



Cervical arthrosis p. 45

Acute sinusitis p. 46

Cerebral vascular disease

Chronic sinusitis p. 47

Classic headache

Facial algia

Migraine

Tonsillitis

Trigeminal neuralgia

Muscle-tensive headache p. 48



Cervical arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	35	20	5



Acute sinusitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	20	10	5

Cerebral vascular disease

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	30	15	5

Chronic sinusitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
15	30	15	5

Classic headache

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	40	10	5

Facial algia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	40	10	5

Migraine

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
8	20	10	5

Tonsillitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
4	20	10	5

Trigeminal neuralgia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
50	40	15	5



Muscle-tensive headache

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	30	15	4 - 5

body



body

body

Dorsal arthrosis	p. 51
Scheuermann morbus	
Humerus-scapular arthrosis	p. 52
Humerus-scapular periarthritis	
Brachial neuralgia	p. 53
Bronchial asthma	
Asthenia	p. 55
Chronic myocardial-hepatopathy	
Insomnia	
Metabolism diseases	
Neuro-vegetative dystonia	
Backache	p. 57
Lumbosacral arthrosis	
Nephrosclerosis	
Colitis	p. 59
Crohn morbus	
Dyspepsia	
Gastralgia	
Gastritis	
Gastroduodenal ulcer	
Ankylosis spondylitis	p. 61
Hip arthrosis	
Sciatalgia - Ischialgia	
Intercostal neuralgia	p. 62



Dorsal arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
25	40	20	4

Scheuermann morbus

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
50	80	20	4



Humerus-scapular arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
25	30	15	4

Humerus-scapular periarthritits

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	35	15	4



Brachial neuralgia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	20	4

Bronchial asthma

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
4	20	15	4



Asthenia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	65	20	4

Chronic myocardial-hepatopathy

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	30	15	4

Insomnia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	35	20	4 - 5

Metabolism diseases

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
15	60	20	4

Neuro-vegetative dystonia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	30	15	4 - 5



Backache

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	60	15	3

Lumbosacral arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	15	3

Nephrosclerosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	80	20	3



Colitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	30	20	3

Crohn morbus

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	35	15	3

Dyspepsia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	50	20	3

Gastralgia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	30	10	3

Gastritis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	20	20	3

Gastroduodenal ulcer

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
4	15	25	3



Ankylosis spondylitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	50	15	3

Hip arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	30	3

Sciatalgia - Ischialgia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	50	20	3



Intercostal neuralgia

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	50	15	4

limbs



limbs


limbs

Epicondylitis p. 65

Peripheral circulatory diseases p. 66

Knee arthrosis p. 67

Hand/foot arthrosis p. 68

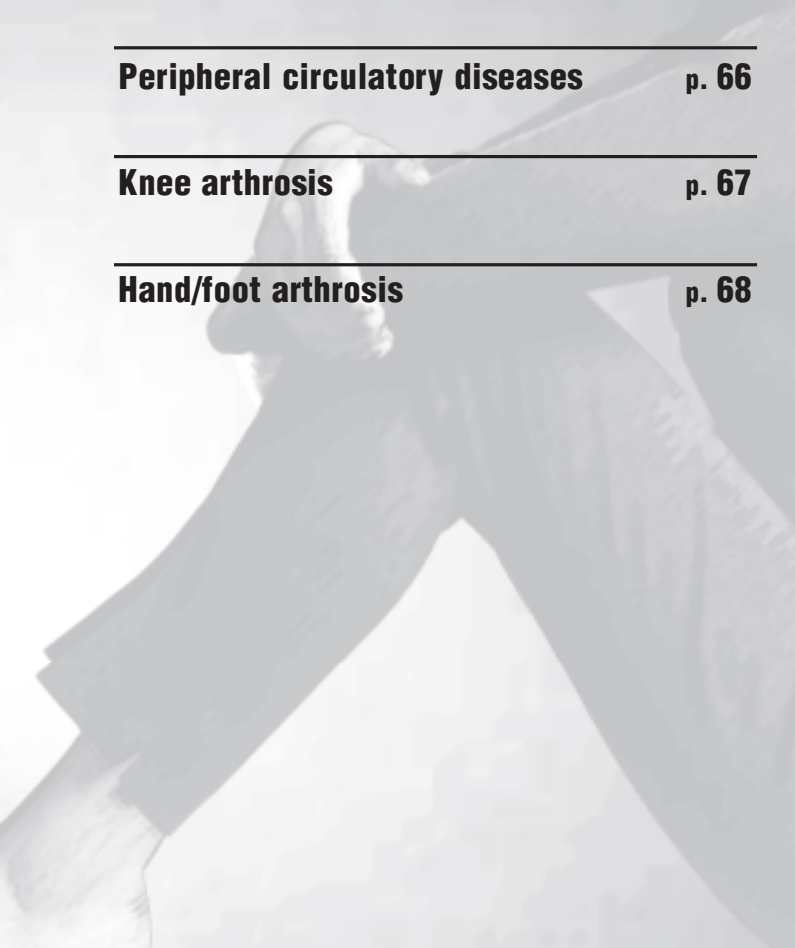



Epicondylitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	30	15	4



Peripheral circulatory diseases

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
15	25	20	2



Knee arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	50	20	2



Hand/foot arthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	35	20	3 - 1

Notes: place solenoid in position 3 for hands and position 1 for feet.

general



general

general

Bursitis p. 71

Burns

Contact dermatitis

Contusions

Convalescence

Fractures

Haematoma p. 72

Hypertonus

Luxation

Muscle tear

Muscle strains

Myositis

Old - age diabetes p. 73

Osteoporosis

Pareses

Phantom limb pain

Pseudoarthrosis

Rheumatoid arthritis

Sprains p. 74

Tendinitis

Ulcus cruris

Wounds

Bursitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	20	15	I.A.

Burns

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
50	60	20	I.A.

Contact dermatitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
8	40	15	I.A.

Contusions

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	50	15	I.A.

Convalescence

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	70	15	I.A.

Fractures

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	50	20	I.A.

I.A.: position solenoid and/or applicators in the Interested Area

Haematoma

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
15	40	15	I.A.

Hypertonus

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	30	10	I.A.

Luxation

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	15	I.A.

Muscle tear

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	50	15	I.A.

Muscle strains

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	60	15	I.A.

Myositis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
5	40	25	I.A.

I.A.: position solenoid and/or applicators in the Interested Area

Old - age diabetes

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	50	20	I.A.

Osteoporosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	30	I.A.

Pareses

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	60	20	I.A.

Phantom limb pain

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	20	I.A.

Pseudoarthrosis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
50	70	20	I.A.

Rheumatoid arthritis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
20	40	15	I.A.

I.A.: position solenoid and/or applicators in the Interested Area

Sprains

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	50	20	I.A.

Tendinitis

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
7	25	15	I.A.

Ulcus cruris

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
10	60	20	I.A.

Wounds

Freq. (Hz)	Int. (%)	Time (min)	Solenoid position
33	60	20	I.A.