



» MEDICART DRAWERS

The Medicart medication trolley with drawers has a sleek and unique look. It has been manufactured from carefully selected hard-wearing materials that have passed the relevant quality tests.

It is a multifunctional, modular trolley with numerous possibilities in terms of structure and applications. The trolleys are amongst other used in A&E departments, by nursing staff, anaesthetists, in theatres...

Our Medicarts have been designed in accordance with the most stringent safety standards and meet all the ergonomic and hygiene requirements.

The trolley's height, the various drawers and the accessories can be tailored to suit your individual needs. It is compatible with standard Euronorm baskets and trays.



» RUNNERS

The runners have been fitted with a **soft-closing system**. The drawer slows down as it closes before being pulled to a complete close in slow motion.


The runners are encased by the drawer frame which makes them invisible and adds to the trolley's sleek design. Contact with grime or grease are a thing of the past. As the drawer is **105% extendible**, trays can easily be lifted out of the frame.



» HANDLES

The drawers have been fitted with aluminium handles, which, thanks to their rounded corners, can be easily cleaned and disinfected.

The handles have been coated with a lacquer that feels soft to the touch.

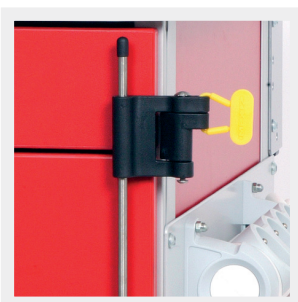

» LOCKS

The Medicart can also be fitted with a lock if required. We propose manual blocking by means of breakable seal.



» KYDEX WORKTOP

The sleek, continuous Kydex worktop comes with raised edges on three sides. As it is devoid of any sharp corners the worktop is easy to clean.

The worktop is scratch proof, can be cleaned with ordinary cleaning products and disinfectants, and is UV and heat-resistant. The surface can be cleaned to perfection, even after coloured disinfectants, e.g. ISO-Betadine, have been spilled.



» CORNER PROFILES

The corner profiles have been designed in such a way that all the Belintra accessories can be fitted without any problems, even by your own technical staff.



» PANELS

The panels are extremely durable, are scratch proof, easy to clean and have a high strength-to-weight ratio.

Thanks to its closed cell structure, the material is ideal for the medical sector.



» CASTORS

The castors meet the highest quality standards and are the product of years of expertise. Studies have shown that double castors reduce rolling resistance by 50% and make the cart easier to manoeuvre.

The Medicart comes with double Tente swivel castors, 1 of which has been fitted with a brake and 1 with a swivel lock (antistatic).

» DOUBLE CASTORS STUDY: LOW ROLLING RESISTANCE AND EASY TO MANOEUVRE

Type of medical cart		Competitor cart 1	Competitor cart 2	Belintra Test cart
Start forward thrust Castors straight	Measurement 1	5,38	7,90	2,79
	Measurement 2	6,40	6,30	3,98
	Measurement 3	7,40	6,90	3,92
	Average	6,39	7,03	3,56
Start forward thrust Castors transverse	Measurement 1	8,30	10,20	7,70
	Measurement 2	8,10	9,50	7,40
	Measurement 3	9,80	12,40	6,10
	Average	8,73	10,70	7,07

MEDR/PS/EN/0001



» MEDICART ROLLER SHUTTER 1100

The Medicart medication trolley with roller shutter is, just like the Medicart with drawers, a multifunctional, modular cart, which distinguishes itself in terms of choice of materials, ergonomics, hygiene, durability and its unique design.

The cart can be fitted with numerous internal components: from baskets to trays, in grid panels or on telescopic runners. They have been designed with your specific needs in mind.

Whether used for medical care or to transport and/or distribute medication, the Medicart is the most ergonomic, hygienic and patient-safe application.

The cart's contents and accessories can be fully tailored to your individual needs and you also have the option of combining baskets and medication trays. It is compatible with standard Euronorm baskets and trays. Special configurations require +12 weeks lead time for delivery.



» PLASTIC ROLLER SHUTTER

The cart is fitted with a push button to ensure that it opens up easily and ergonomically. As the roller shutter retracts into a sealed unit in one go, dirt can never play havoc with the system.

Thanks to the design of the cart, both the inside and the outside of the roller shutter are easy to clean.



» HANDLE

The plastic roller shutter is fitted with an aluminium grip that runs across its entire width. This grip ensures that the roller shutter can be opened and closed again safely and ergonomically.



» LOCK

The shutter comes with an electronic number lock which blocks the push button when closed and can be opened by means of a programmable 4-digit code. Both the push button and the electronic number lock are located directly underneath the worktop. One set of batteries will allow you to operate the electronic number lock 80,000 times. If the batteries are dead, the cart can always be opened with an emergency key.



» KYDEX WORKTOP

The sleek, continuous Kydex worktop comes with raised edges on three sides. As it is devoid of any sharp corners the worktop is easy to clean.

The worktop is scratch proof, can be cleaned with ordinary cleaning products and disinfectants, and is UV and heat-resistant. The surface can be cleaned to perfection, even after coloured disinfectants, e.g. ISO-Betadine, have been spilled.



» CORNER PROFILES

The corner profiles have been designed in such a way that all the Belintra accessories can be fitted without any problems, even by your own technical staff.



» PANELS

The panels are extremely durable, are scratch proof, easy to clean and have a high strength-to-weight ratio.

Thanks to its closed cell structure, the material is ideal for the medical sector.



» CASTORS

The castors meet the highest quality standards and are the product of years of expertise. Studies have shown that double castors reduce rolling resistance by 50% and make the cart easier to manoeuvre.

The Medicart comes with double Tente swivel castors, 1 of which has been fitted with a brake and 1 with a swivel lock (antistatic).

» DOUBLE CASTORS STUDY: LOW ROLLING RESISTANCE AND EASY TO MANOEUVRE

Type of medical cart		Competitor cart 1	Competitor cart 2	Belintra Test cart
Start forward thrust Castors straight	Measurement 1	5,38	7,90	2,79
	Measurement 2	6,40	6,30	3,98
	Measurement 3	7,40	6,90	3,92
	Average	6,39	7,03	3,56
Start forward thrust Castors transverse	Measurement 1	8,30	10,20	7,70
	Measurement 2	8,10	9,50	7,40
	Measurement 3	9,80	12,40	6,10
	Average	8,73	10,70	7,07

MER511/PS/EN/0001

» ACCESSORIES FOR THE MEDICART DRAWERS

With a myriad of accessories to choose from, the Medicart can be equipped to fit the needs of every department.

Please note that not all accessories fit every position on the Medicart and that the number of accessories which can be mounted on a Medicart is limited.



» GREY METAL BIN

Status	Article Number	Description	Capacity	Ext. Dimensions
A	19010 00006	-	18 l	300 x 175 x 360 mm
B	19010 00017	2 next to each other	2 x 18 l	440 x 255 x 380 mm



» WALLY® WASTE BIN WITH KNEE LEVER

Status	Article Number	Description	Capacity	Ext. Dimensions
A	19010 00023	-	8 l	230 x 210 x 375 mm
A	19010 00024	2 next to each other	2 x 8 l	465 x 210 x 375 mm



» FOLD-DOWN WORKTOP

Status	Article Number	Maximum weight	Ext. Dimensions
A	19010 00008	25 kg	505 x 400 x 15 mm



» DRAIN FIXTURE

Status	Article Number	Description	Int. Dimension	Ext. Dimensions
A	19010 00001	For Stala structure	132 mm	600 x 134 x 82 mm
B	19010 00002	Side drain for Medicart	132 mm	500 x 137 x 85 mm
A	19010 00014	Side drain for Medicart	160 mm	500 x 165 x 85 mm

» MEDICART DRAWERS ACCESSORIES



» STALA BASE STRUCTURE FOR DIFFERENT TRAY TYPES

Status	Article Number	Description	Ext. Dimensions
B	19010 00010	For Medicart drawers	Structure without trays 652 x 20 x 650 mm



» TRAYS FOR STALA STRUCTURE

Status	Article Number	Description	Status	Article Number	Description
B	15010 00001	3 trays	A	15010 00002	4 trays
B	15010 00004	6 trays	A	15010 00003	5 trays
B	15010 00005	9 trays			



» FIFO PICBOX FOR STALA STRUCTURE

Status	Article Number	Description	Ext. Dimensions
B	15010 00008	3 trays	600 x 200 x 245 mm
B	15010 00009	5 trays	600 x 135 x 165 mm



» DIN-RAIL

Status	Article Number	Description	Ext. Dimensions
A	19010 00009	For the side of Medicart	510 x 25 x 10 mm
B	19010 00029	For Stala structure	620 x 25 x 10 mm



» CATHETER HOLDER FOR DIN-RAIL

Status	Article Number	Ext. Dimensions
A	30004 00001	135 x 85 x 530 mm

» MEDICART DRAWERS ACCESSORIES



» CONTAINERHOLDER FOR NEEDLES

Status	Article Number	Ext. Dimensions
A	19010 00007	120 x 121 x 200 mm



» CONNECTING UNIT FOR PC BASKETS OF 300 X 400 MM

Status	Article Number	Ext. Dimensions
A	19010 00003	500 x 345 x 140 mm



» FLEX ARM FOR PORTABLE

Status	Article Number	Description	Maximum weight	Ext. Dimensions
B	14009 00001	-	8 kg	370 x 290 mm
B	14009 00037	With clamps for portable	8 kg	420 x 250 x 30 mm



» KENSINGTON LOCK FOR PORTABLE MEDIA

Status	Article Number	Description
B	14009 00013	With number lock



» CPR BOARD


Status	Article Number	Ext. Dimensions
A	19010 00031	540 x 450 x 15 mm

» MEDICART DRAWERS ACCESSORIES



» OXYGEN BOTTLE HOLDER

Status	Article Number	Diameter
A	19010 00045	Universal dimension

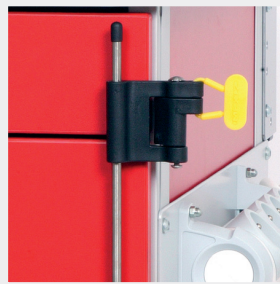

» SEAL LOCK

Status	Article Number	Description
A	19010 00034	Lock
A	19010 00041	Seals (100/pack)



» ELECTRIC SOCKETS WITH EARTH PINS CEE 7/5 TYPE E

Status	Article Number
A	19010 00036



» STAINLESS STEEL SERUM HOLDER

Status	Article Number	Description
A	14008 00004	With 2 hooks
B	14008 00006	With 4 hooks
A	14008 00011	Corner unit in angle of 90°



» DEFIBRILLATOR HOLDER

Status	Article Number	Ext. Dimensions
A	19010 00033	370 x 400 mm

» MEDICART DRAWERS ACCESSORIES



» CONNECTING UNIT FOR 3 GLOVE BOXES

Status	Article Number	Int. Dimensions	Ext. Dimensions
A	19010 00039	90 x 252 x 418 mm	255 x 93 x 420 mm