



# dicoNEX SN Single-use syringe, 3-part, luer, with needle alongside the syringe, sterile

att.no.3 pt.2 subpt.2.3

revision 2

date: 2019-11-12

TD-41

amendment of: 2023-02-01

**MANUFACTURER:** ZARYS International Group sp. z o.o. sp.k.


**ADDRESS:** ul. Pod Borem 18, 41-808 Zabrze, Poland, phone +48 32 376 07 00-02, fax. +48 32 370 38 08

## INDICATION:

- 3-piece syringe for subcutaneous, intramuscular and intravenous infusions
- Intended for drawing up and administering of medicines / liquids to patient and collecting body fluid samples from patient

## FEATURES:

- in set with needle alongside the syringe
- transparent barrel allows full visual control of content
- green contrasting plunger
- luer nozzle placed centrally for 2 ml capacity, eccentric for 5, 10, 20, 50 ml capacity
- sealing in the form of double ring on plunger stopper situated on the top of plunger
- smooth plunger sliding thanks to rubber sealing
- limitation flange securing against accidental removal of plunger
- name of syringe and manufacturer's logo placed on barrel
- clear, legible and permanent black graduation facilitating dosage, adequately:
  - every 0,1 ml for capacity of 2 ml
  - every 0,2 ml for capacity of 5 ml (extended to 6 ml)
  - every 0,2 ml for capacity of 10 ml
  - every 0,5 ml for capacity of 20 ml
  - every 1 ml for capacity of 50 ml
- needle tube is made of smooth and durable stainless steel
- needle tip sharpened in three surfaces
- bevel type: lb/bl standard long bevel
- hubs coloured in accordance with iso code
- single-use, non-toxic, non-pyrogenic
- latex-free
- phthalate-free
- EO sterilized
- shelf life: 5 years
- packaging: 1 pc./blister-pack; 25/50/100 pcs/box

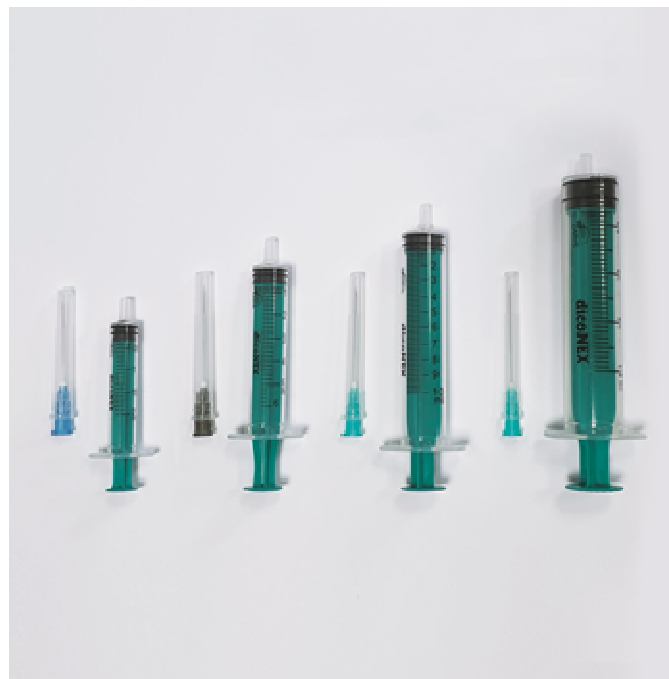

## SIZE AND REFERENCE

CATALOGUE NUMBER	CAPACITY	NEEDLE SIZE	NEEDLE COLOUR	INDIVIDUAL PACKAGING	INTERMEDIATE PACKAGING	BULK PACKAGING
002ML-3CZ-SN	2 ml	0,6mm x 30mm	deep blue	blister-pack	100 pcs	30 x 100 pcs
005ML-3CZ-SN	5 ml	0,7mm x 30mm	black	blister-pack	100 pcs	24 x 100 pcs
010ML-3CZ-SN	10 ml	0,8mm x 40mm	deep green	blister-pack	100 pcs	16 x 100 pcs
020ML-3CZ-SN	20 ml	0,8mm x 40mm	deep green	blister-pack	50 pcs	16 x 50 pcs
050ML-3CZ-SN	50 ml	0,8mm x 40mm	deep green	blister-pack	25 pcs	16 x 25 pcs